

AiveoHP

Balloon Dilatation Catheter

The world's smallest profile
PCI Balloon Catheter



Technical Specification	
Balloon Diameter (mm)	0.75 、 1.0 、 1.25 、 1.5 、 1.75 、 2.0
Balloon Length (mm)	5 、 8 、 10 、 12 、 15 、 20 、 25 、 30
Nominal Pressure	10 atm
Rated Burst Pressure	20 atm
Average Burst Pressure	≥ 28 atm
Compliance	Semi-compliant
Balloon Pleats	Bi-fold (2 pleats) for Ø0.75-1.0mm and Tri-fold (3 pleats) for Ø1.25-2.0mm
Tip Length	1.5mm for Ø0.75-1.75mm and 2.0mm for Ø2.0mm
Proximal Shaft	1.9F (Ø0.75-1.5mm) / 2.0F (Ø1.75-2.0mm)
Distal Shaft	2.2F (Ø0.75-1.0mm) 2.2F/2.36F (Ø1.25-1.5mm) 2.36F/2.55F (Ø1.75-2.0mm)
Working Length (cm)	140
Guide Wire Compatibility	0.014 inch
Guiding Catheter Compatibility	5F Compatible
Coating	•Slider™ Hydrophilic external surface coating from tip to adjacent guidewire port •Eel™ Hydrophobic coating inside guidewire lumen

Ordering Information								
Balloon Diameter	Balloon Length							
	5mm	8mm	10mm	12mm	15mm	20mm	25mm	30mm
Ø 0.75mm	813-0705	813-0708	813-0710	813-0712	813-0715	813-0720	813-0725	813-0730
Ø 1.0mm	813-1005	813-1008	813-1010	813-1012	813-1015	813-1020	813-1025	813-1030
Ø 1.25mm	813-1205	813-1208	813-1210	813-1212	813-1215	813-1220	813-1225	813-1230
Ø 1.5mm	813-1505	813-1508	813-1510	813-1512	813-1515	813-1520	813-1525	813-1530
Ø 1.75mm	813-1705	813-1708	813-1710	813-1712	813-1715	813-1720	813-1725	813-1730
Ø 2.0mm	813-2005	813-2008	813-2010	813-2012	813-2015	813-2020	813-2025	813-2030

For more information please visit our website at www.brosmed.com or contact us



BrosMed Medical Co., Ltd
15th Building,SMEs Venture Park SongShan
Lake Hi-Tech Industrial Development Zone
Dongguan 523808 ,GuangDong ,China
Tel: +86-769-2289 2018
Fax: +86-769-2289 2016
E-mail: sales@brosmed.com

BrosMed Medical B.V.
Mgr. Buckxstraat 8, 6134 AP Sittard ,The Netherlands
Tel: + 31 850 140 901
Mail: CS@BrosMed.com
Web: www.brosmed.com

CE 2292



AiveoHP

Balloon Dilatation Catheter



CTO HP Ø 0.75

- Smallest Balloon Diameter
- Smallest Crossing Profile
- High Pressure Dilatation

www.brosmed.com


Alveo HP

Balloon Dilatation Catheter



Ultra-low Tip Profile

Leading laser processing technology ensures a smooth, tapered ultra-low tip profile of **0.0156"**

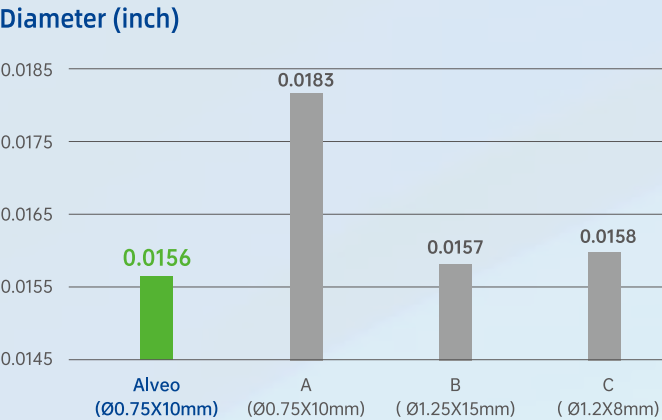

Entry Profile
0.0156"

* Competitor data on file.




Designed for CTO and Complex Lesions

RBP: **20atm** ABP: **≥ 28atm**

* In Vitro lab data. Do not exceed RBP.


Balloon Diameter: **0.75mm**
Rated Burst Pressure: **20atm**

Robust balloon material allows high pressure and improved kink resistance for CTO and complex lesions.

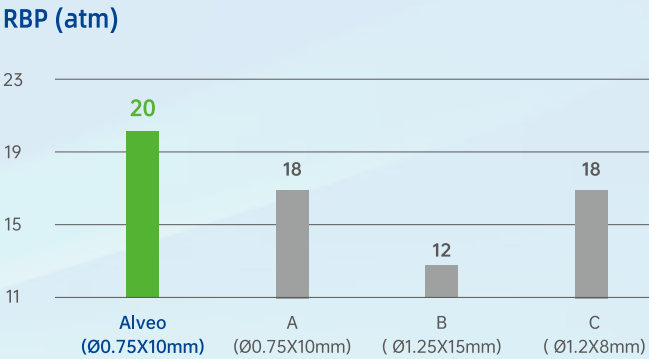

* Competitor data on file.

Balloon Profile

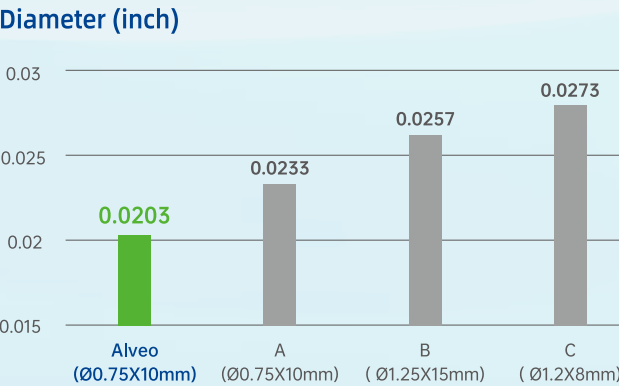

Note : marker band OD **0.020"**
* Competitor data on file.

Folded Balloon Profile
0.0203"



0.0184"
Ultra-low Crossing Profile